Yoghi

MASS BODY 550

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EN

## Safety Instructions

BEFORE USING THE DEVICE, PLEASE FAMILIARIZE YOURSELF WITH ALL THE SAFETY INSTRUCTIONS AND INSTRUCTIONS FOR USE! PASS ON ALL DOCUMENTS CONCERNING THE DEVICE WHEN YOU GIVE IT TO A THIRD PARTY!

In the event of damage resulting from non-compliance with the operating instructions, the warranty claim is void! No liability is assumed for consequential damage! No liability is assumed for damage to property and persons resulting from improper use or failure to observe the safety instructions!


- **This device may be used by children over 8 years of age and persons with limited physical, sensory or mental abilities, or lacking experience and knowledge, provided they receive supervision or instructions regarding safe use. device safety and understand the risks it presents. It is forbidden to let children play with the device. It is prohibited to allow children to clean and maintain the appliance without supervision.**
- **Keep the unit and its cable out of the reach of children.**
- **Never leave children unattended with packing materials. Packaging materials are a suffocation hazard. Children often underestimate the resulting dangers. Always keep children away from packaging materials.**
- **Do not repair the unit yourself. Contact a qualified technician. To avoid danger, has the defective cable replaced only by an equivalent cable and only by the manufacturer, its service department or a similarly qualified person.**
- **Always disconnect the unit from the power supply if left unattended and before assembly, disassembly or cleaning.**
- **Turn off the unit and disconnect it from the power supply before changing accessories or approaching parts that move during operation.**

- This unit is not designed to be operated by an external timer or separate remote-control system.
- **WARNING!** Always keep the unit dry.
- Never touch the power plug with wet or damp hands.
- Do not operate it in damp or humid conditions.
- Never immerse the unit in water or any other liquid. Do not run the parts under running water.
- Unplug the AC adapter cable from the wall outlet: when a malfunction occurs, before cleaning the product, during a thunderstorm.
- When unplugging, remove the power plug directly from the outlet but do not pull on the power cord.
- Never insert pins or other metal fasteners into the unit.
- Do not use the unit if it shows visible signs of damage. Do not use any damaged unit.
- If liquid enters the unit, check it before using it again.
- Do not continue to use the unit if you notice any damage to the power cord, plug, housing or other parts.
- Keep the AC adapter cable away from heated surfaces.
- Excessive use may cause overheating and shorten the life of the unit.
- Do not operate under a blanket or pillow. Doing so may cause risk of fire, electric shock or injury.
- Make sure that the massage mechanism or any other moving part of the product is not obstructed by hair, clothing or jewelry during operation.
- Repairs should be made only by qualified service personnel. Improper repairs may result in considerable risk to the user.
- Defective parts should always be replaced with original spare parts.
- Position the AC adapter cable so that it cannot be stepped on or tripped over.
- Do not use the device for head massage.
- Do not use it on animals.

- Please do not use it on people with herniated discs or abnormal spinal changes.
- Never use the device to massage the heart area of people who are equipped with it. Always consult a physician before using the device on other parts of the body of a person with a pacemaker.
- Never use the device on swollen, burned, inflamed or otherwise injured skin or body parts.
- The device should not be used during pregnancy.
- Do not use the device on the face (eyes), larynx or other particularly sensitive parts of the body.
- Never use the device while you are sleeping.
- Never use the device after taking medication or after consuming alcohol (reduced perception).
- This is a non-professional device, intended for personal use and designed to provide a soothing massage to tired muscles. Do not use as a substitute for medical treatment.
- Before using the massager, you should consult your physician if you have a serious illness or are recovering from surgery.
- Protect the unit from moisture, dripping and splashing.
- Do not place the unit on or near hot surfaces or near flames.
- Do not use any corrosive or scouring cleaning agents.
- **CAUTION:** Removal/disposal of rechargeable batteries must be done by the manufacturer, its service personnel or similarly qualified persons to avoid a hazard.
- The battery must be disposed of safely.
- **WARNING:** Unplug before removing the battery!
- Take the battery to a collection point for batteries, cells and/or hazardous waste.

## Parts description

- |                                  |                                  |
|----------------------------------|----------------------------------|
| 1. LED (operation indicator)     | 8. Acupuncture massage head      |
| 2. Key (mode selection)          | 9. Convex massage head           |
| 3. Key (increase power level)    | 10. The mowing massage head      |
| 4. Button (decrease power level) | 11. The four-legged massage head |
| 5. Key (on/off)                  | 12. The anti-dust latch          |
| 6. Handle                        | 13. The charging port            |
| 7. Thread                        | 14. AC adapter cable             |
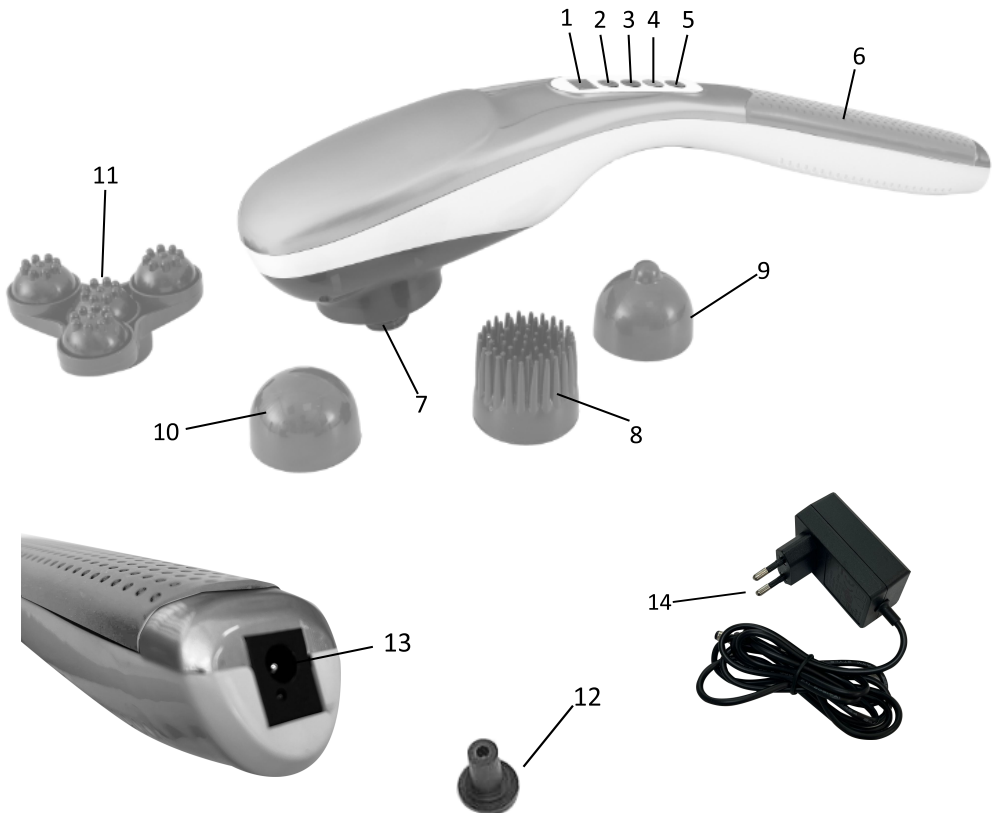
## Technical data

Adapter voltage: 100-240V~, 50/60Hz

Power consumption: 16 W

Speed setting of 3000-3200 rpm

Rechargeable lithium battery: 1800 mAh

Main material: ABS

## Before using for the first time

- All parts must be unpacked and the packaging materials must be completely removed.
- Make sure that the contents of the package are complete and intact!

## Use

The dust latch must be removed to use or charge the product.

Charging: Remove the dust cover.

Connect the AC adapter cable to the unit.

Plug the adapter into a wall outlet.

The indicator light will turn red and will turn green when the unit is fully charged.

Removing the massage head: Turn one massage head counterclockwise.

Attach the massage head: Screw on a massage head by turning clockwise.

Turn on the device: Press . A signal sounds. The red LED is permanently lit. The device is in mode 1 on power level 1.

Switching off the device: Press . A signal sounds. The red LED goes out.

Settings:  Switch between modes to set different massage styles. 6 modes are available.

 Increase power level

 Decrease power level

## **Cleaning and Maintenance**

**DANGER!** The risk of electric shock!

Before cleaning the product, disconnect the power plug from the wall socket.

Never touch the adapter with wet or damp hands. Always keep the adapter and power cord dry. Turn the unit off.

Never immerse the product in water or any other liquid. Never hold the unit under running water.

Do not use any corrosive or scouring cleaning products.

Clean the unit and massage heads with a damp cloth. For stubborn dirt, apply a few drops of mild detergent to a damp cloth and clean the unit again.

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