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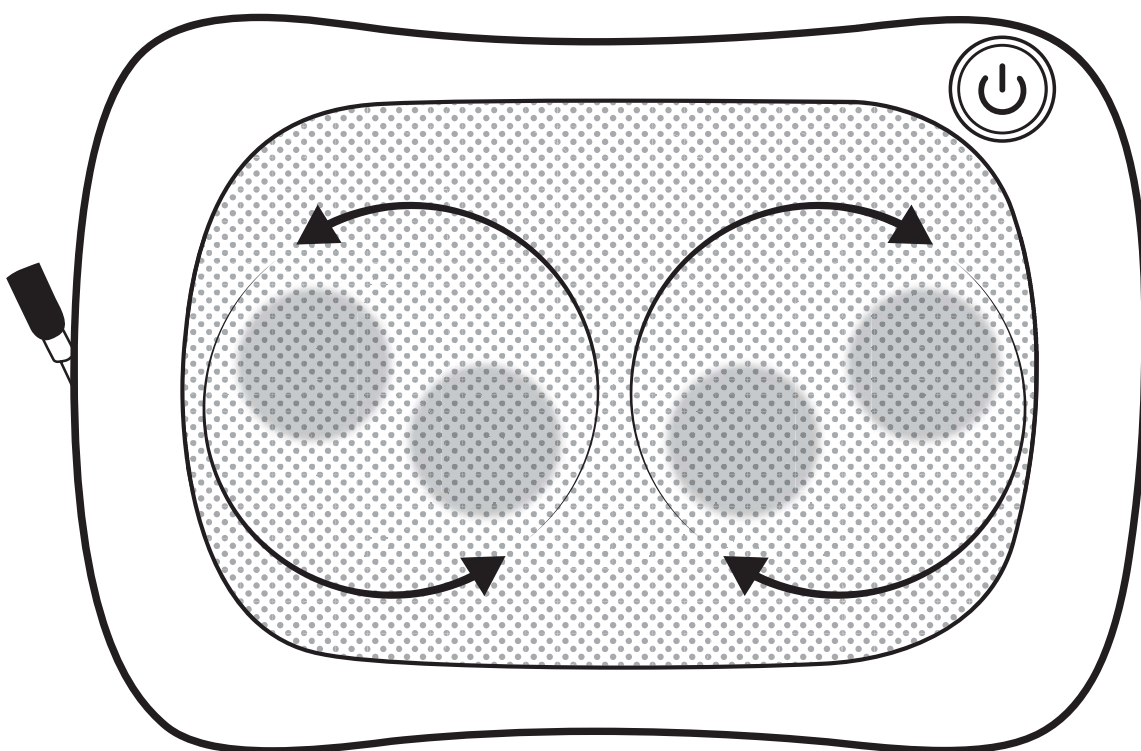
Dear Customer,

Thank you for purchasing the Multi Massager. Please read the instructions before using the massager as this manual provides additional information including safety precautions and warnings.

If you have any difficulties or need more information, you can find our help centre by scanning the QR code opposite with your smartphone. If you cannot find the answer to your question, you can also contact our After Sales Service: support.uk@terraillon.com

DISCLAIMER

CAUTION

If you have the following conditions or any health issues, please consult with your doctor before use.

1. Pregnancy, diabetes, osteoporosis.
2. Use of cardiac pacemaker, recent surgery, epilepsy or migraines, protrusion of lumbar intervertebral discs and anterior displacement of bone vertebrae.
3. Metal implants and other concerns of the body.
4. In vivo transplanted medical electronic devices such as pulse regulators.
5. Artificial heart, lungs and other life-sustaining medical electronic devices.
6. Wearable medical electronic devices such as electrocardiographs.
7. Patients with acute disease.
8. Patients who feel abnormalities in the body.

PRECAUTIONS


1. Minors are not allowed to use this product. This product is not a toy, please store it out of reach of children.
2. The device has a heating function, and people who are not sensitive to heat must pay attention to it. Patients with infectious diseases should use this product alone. If you notice discomfort, discontinue use immediately.
3. This product is a personal care appliance and cannot replace a medical device.
4. Do not immerse the product in water. If the product falls into the water, please submit it for inspection in time.
5. Do not place the product near a heat source or into the fire.
6. Do not drop the product intentionally.
7. Do not attempt to disassemble the product.
8. It is forbidden to replace or supplement medical treatment with the device, if you feel pain or are receiving treatment, please consult before using the device your doctor.
9. If the equipment is damaged, stop operation immediately.
10. Do not use in places such as bathroom.
11. Do not use this product while intoxicated or unwell.
12. Do not fall asleep while the product is being handled.
13. It's forbidden to use product accessories from other suppliers.

PACKAGE CONTENTS

Multi Massager x 1
Instructions Manual x 1
Charging cable x 1

PRODUCT OVERVIEW

Shiatsu Multi functional massage pillow

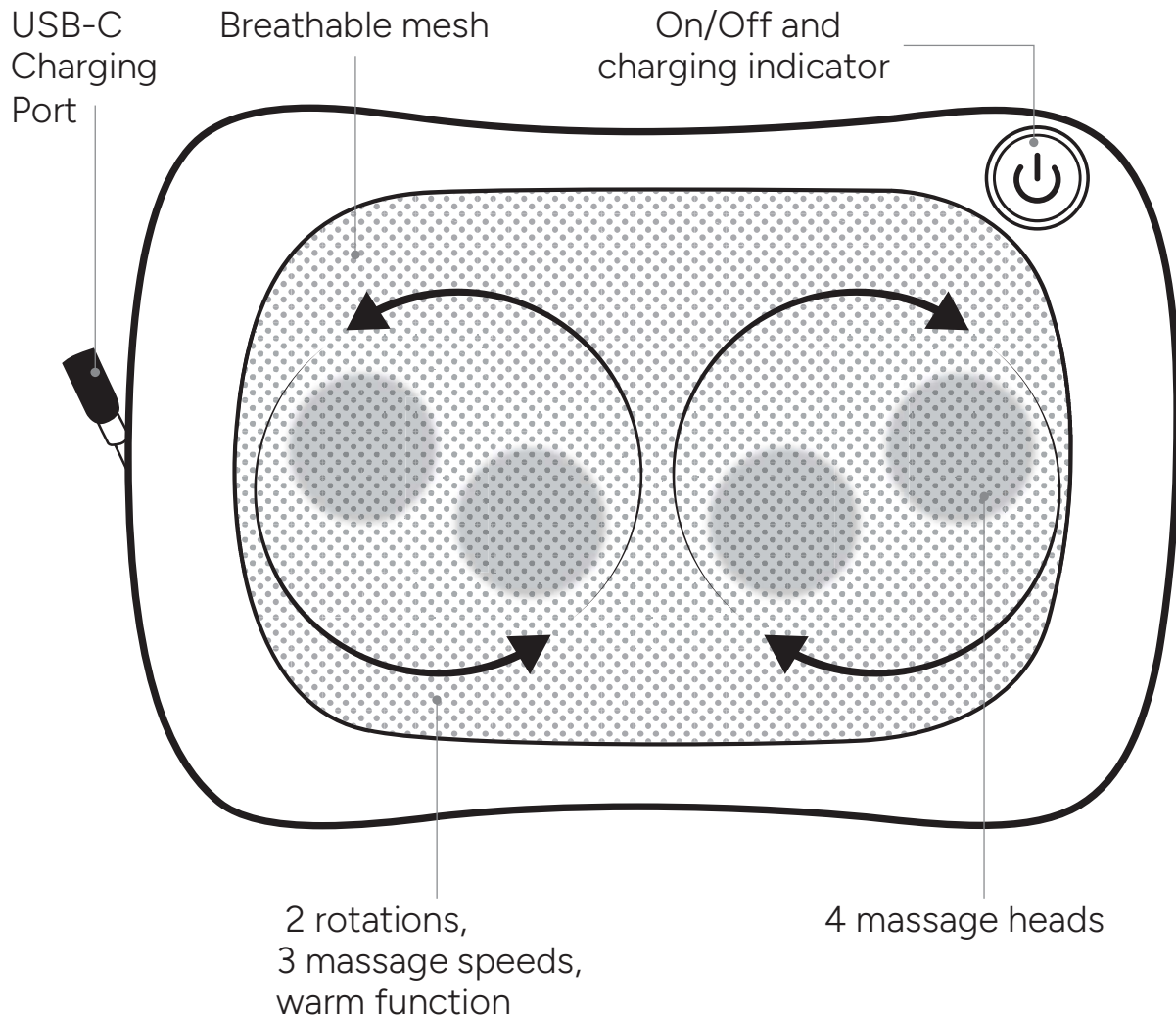CHARGING


For your safety, always use an appropriate 7V input charger/adaptor to charge the massager. DO NOT USE a quick-charging USB adapter as this may damage the product.

Important: the Multi Massager must not be turned on while charging.

1. Use the USB-C cable included with this massager.
2. Plug the USB-C plug into the USB-C charging port.
3. Connect the USB plug to a USB port on a standard wall adapter to start charging.
4. When the Multi Massager is charging, the indicator light is red. When fully charged, the indicator light will turn green.

NOTE:

- You must select a qualified charger from a reputable manufacturer.

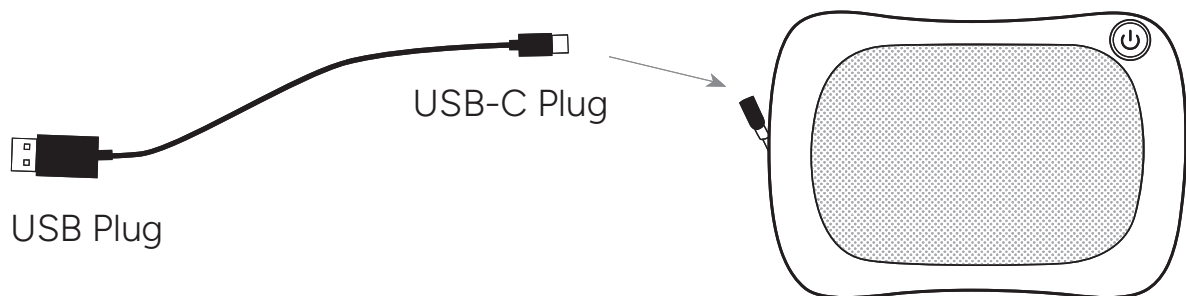 Safety instructions for charging the product's lithium battery:

- The input voltage is DC 7V 1A.
- In order to charge the product safely and reliably, please use a wall adapter that is certified and in good condition.
- Only charge the product in ambient temperatures between 10°C /40°F and 40°C /104°F.
- Do not heat above 70°C /158°F. Do not dispose of in a fire.
- Keep out of direct sunlight.
- Turn off the product after use.

- Please note we will not take any responsibility for the misuse of the product, such as using voltage higher than indicated. Misuse may result in serious injury or loss of property.
- In extreme cases, abuse or misuse of a lithium battery can lead to an explosion, heat generation, fire or smoke development.

USING THE PRODUCT

1. Press the on/off button  once to turn on the device at low speed and turn heating forward. The heating function will be activated by a red light and the massage heads will begin to rotate.
2. Press twice to activate the medium speed.
3. Press three times to activate the high speed.
4. Press four times to turn off.

During the massage cycle, the direction of rotation of the heads changes every minute.

The duration of a program is approximately 15 minutes. The appliance will stop automatically at the end of the program.

FAQ

• Why is my Multi Massager not turning on?

If the Multi Massager will not power on, please ensure the battery is fully charged.

• Why is my Multi Massager is not charging?

If your Multi Massager will not charge, please ensure each end of the charging cable is securely attached and connected to the power source. Check to ensure the wall outlet has power. If using a computer USB port, please ensure that the computer has power.

TECHNICAL SPECIFICATIONS

Product name	Multi Massager
Battery capacity:	2 000 mAh
Massage function	4 massage heads, 2 rotations, 3 massage speeds
Warm function	Breathable mesh
Charging time	4h30
Fully charged battery usage	About 1h30
Charging port	USB-C
Output power	18 W
Product size	33 cm
Cable size	120 cm
Weight	1,6 kg

CARE AND MAINTENANCE

To reduce the risk of failure, electric shock, injury, fire and equipment damage, please read all the instructions and warnings in this manual.

1. Do not dispose of, store or use the massager near or in fire or open flames.
2. Do not attempt to disassemble the massager.
3. Please consult with your doctor before use. This massager is not intended to replace or supplement prescribed medical treatments.
4. Only wipe clean with a damp towel, don't use corrosive detergents.

WARRANTY

This product is guaranteed against defects in materials and workmanship. This warranty does not cover damage resulting from accidents, improper use or negligence. If you have a complaint, first contact the store where you made your purchase.

ENVIRONMENTAL PROTECTION

When the product is no longer usable, please take it to a collection point for the recycling of waste electrical and electronic equipment.

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Terraillon SAS France & Headquarters
Terraillon Asia Pacific Ltd
Terraillon Corp USA