
Yoghi

LUNE DE SOIE

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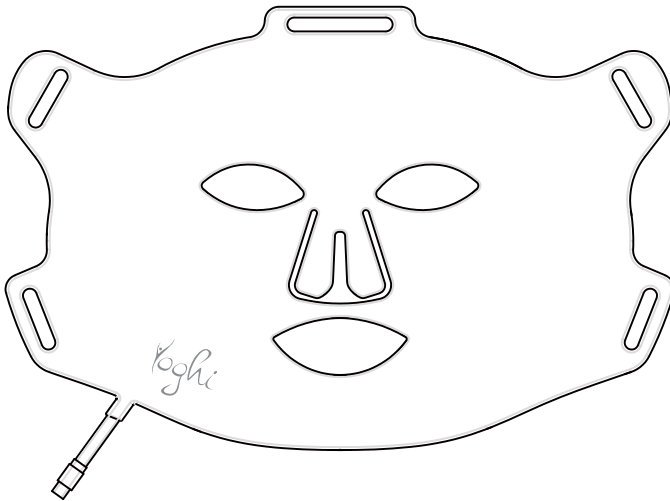EN

SAFETY INSTRUCTIONS

- **This device may be used by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge, if they are properly supervised or have been given instructions concerning safe use of the appliance and understand the risks involved.**
- **This device must only be used by persons aged 18 years or older.**
- **Children must not play with the appliance.**
- **Cleaning and maintenance must not be carried out by children without supervision.**
- **Keep the device and its cable out of reach of children under 8 years old.**
- **Do not allow children to play with packaging materials.**
- **Never attempt to repair the device yourself: contact a qualified technician.**
- **If the power cable is damaged, it must be replaced by the manufacturer or a qualified professional.**
- **This product contains a lithium battery.**
- **Keep the product away from children and pets.**
- **Do not disassemble, modify or repair without authorization.**
- **Never leave the device unattended while in operation.**
- **Do not attempt to switch on the mask while charging: unplug before use.**
- **Do not subject the device to strong impacts or drops.**
- **Do not expose to heat, water, moisture or direct sunlight.**
- **This product is not waterproof: do not immerse in water or rinse.**

- If water or a foreign object enters the device, stop using it immediately.
- Do not use non-original accessories.
- Store the device in a bag or box when not in use, away from dust and sunlight.
- Due to the USB/Type-C cable length, keep it out of reach of children to avoid strangulation risk.
- If you feel discomfort during use, stop immediately.
- Do not use if any part is damaged.
- To avoid water exposure, the LED Light Therapy Mask has an IP22 protection rating.
- Do not use the device in a humid environment (e.g., bathroom, near a shower or pool). Keep it away from water or moisture to prevent electrical hazards or damage.
- When not in use, store the device in its box or carry bag. Keep away from dust, lint and direct sunlight.

PEOPLE WHO MUST NOT USE THIS PRODUCT

Use is strictly prohibited to the following individuals to prevent serious risks:

- Patients with heart diseases.
- Patients with autoimmune diseases.
- Persons with melanoderma.
- Persons with malignant tumors.
- Patients with epilepsy, photosensitive epilepsy or migraines.
- Persons undergoing treatment with photosensitizing medications (photodynamic therapy).
- Persons under 18 years of age.
- Pregnant or breastfeeding women.
- Persons with photosensitivity-related diseases (e.g. lupus

erythematosus, photo-eczema, albinism).

- Persons undergoing dermatological treatment in its initial phase.
- If you suffer from photosensitivity disorders, do not use this product.
- Do not use on the same day as treatments containing photosensitizing ingredients.

PRECAUTIONS BEFORE USE

- If you have a serious illness or are under medical treatment, consult your doctor before use.
- If you have suspicious or cancerous skin lesions, consult a dermatologist before use.
- Skin condition may vary with season, stress or alcohol consumption. If your skin is irritated or abnormal, stop using until recovery or consult a doctor.
- **Skin Sensitivity Test:** Test on a small area before first use. Stop immediately if irritation or allergy occurs.
- **Eye protection:** Keep eyes closed or wear the supplied eye shields. Direct LED light exposure may cause discomfort or eye damage.
- **Clean and dry skin:** Remove makeup and cosmetics before use to avoid interference or irritation.
- **Duration and frequency:** Follow recommended use: 3–4 times per week, sessions of 10–20 minutes.
- This product may not be suitable for certain skin or medical conditions. When in doubt, consult a healthcare professional.

USE WITH CAUTION / MEDICAL CONSULTATION REQUIRED

- Fever or abnormal body temperature.
- People with high blood pressure or blood disorders.

- People with dermatoses or atopic dermatitis.
- Individuals currently under outpatient medical treatment.

DO NOT USE IN THE FOLLOWING AREAS / CONDITIONS

To avoid accidents and severe skin complications:

- Eyeballs.
- Upper eyelids (eye socket).
- Head.
- Areas with pathological pigmentation.
- Viral warts.
- Fracture areas.

For areas recovering from cosmetic surgery, consult a medical professional first.

ADDITIONAL SAFETY

- Do not place any naked flame, such as lit candles, on the device.
- Do not use the device in environments with strong electromagnetic interference.
- **CAUTION:** Do not subject the device to heavy shocks to prevent damage to the lamp.

LIGHT SENSITIVITY TEST

Do not use this product if you fall under the categories listed below:

- Pregnant or planning to become pregnant, or breastfeeding.
- Persons sensitive to light or taking photosensitizing medications; persons with hereditary eye diseases (e.g. lupus erythematosus, photo-eczema, albinism).

Step 1: Place the device on your arm for 10 minutes of irradiation. If redness or itching occurs, you may have photosensitive skin. Stop use or consult a dermatologist.

Step 2: Clean your face and keep it dry before each use (no makeup) to avoid interactions with photosensitive ingredients.

Step 3: Redness may appear after use, which is normal and disappears within 1–2 hours. If burning or discomfort occurs, stop use immediately.

PRODUCT OVERVIEW (FIG.1)

1. LED Light Therapy Mask
2. Remote control
3. Silicone eye shields
4. Storage pouch
5. USB Type-C cable
6. Adjustable straps

LIGHT MODES OVERVIEW

• Yellow light (590 nm)

Balances skin texture, stimulates red blood cells, reduces redness, swelling and fine lines.

• Red light (660 nm)

Anti-aging: boosts collagen production, firms and improves elasticity.

• Near-infrared light (850 nm)

Repairs skin, enhances blood circulation and hydration.

• Blue light (415 nm)

Eliminates acne-causing bacteria, reduces inflammation and scarring, boosts oxygenation.

SIX MODES

- **M1: Repair mode** – 660 nm (red)
- **M2: Rejuvenation mode** – 590 + 660 nm (yellow + red)
- **M3: Anti-aging mode** – 660 + 850 nm (red + near-infrared)
- **M4: Morning care mode** – 590 + 660 + 850 nm
- **M5: Anti-acne mode** – 415 nm (blue)
- **M6: Evening care mode** – 590 + 660 + 850 + 415 nm

BEFORE FIRST USE

Fully charge the mask battery (approx. 3 hours).

Remove the protective film.

Wear the eye shields.

Attach the strap around your head as shown in FIG.2.

Connect the mask to the remote control.

REMOTE CONTROL (FIG.3)

- a. Power port (input/output)
- b. OLED display
- c. Power / mode selection button
- d. Brightness adjustment
- e. Pulse button

OPERATION

Power on: Press and hold the power button for approx. 3 seconds. The light turns on.

Brightness adjustment: Short press to switch brightness level: 50%, 75% or 100%.

Timer: Each mode is preset to 10 minutes.

Mode selection: Press the power button again to cycle through the 6 modes.

Pulse function: Can be activated at any time in any mode.

Charging: Connect the USB Type-C cable and plug into a power source. The battery icon flashes during charging and shows full capacity when complete (approx. 3 h).

After the session

- Apply your usual skincare products (serum, moisturizer).
- LED light therapy enhances product absorption.

Additional precautions

- If concerns or adverse effects occur, consult a healthcare professional or contact customer service.
- This device is not a substitute for medical advice or dermatological care.
- Maximum light surface temperature: 42 °C (safe for 10-minute sessions).

CLEANING, MAINTENANCE AND STORAGE

- Wipe gently with a soft slightly damp cloth or non-abrasive wipe.
- Do not use harsh chemicals or alcohol.
- Ensure the mask is completely dry before storage.
- Store in a clean, dry place away from sunlight and humidity.
- Do not fold the mask or place heavy objects on it.
- Avoid extreme temperatures and humid environments.
- If stored in extreme conditions, allow at least 2 hours before use.

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